

beurer quick start guide

WIRELESS
CONNECT
BG900



Supported browser up to release date:

- Windows® Internet Explorer® 9.0 → Mozilla® Firefox® 11.0
- Safari 5 (Mac) → Google Chrome™ 18.0 → Opera™ 11.0

www.beurer-connect.com



Beurer GmbH
Söflinger Straße 218
89077 Ulm · Germany
Tel. +49(0)731/3989-144
Fax +49(0)731/3989-255
www.beurer.de · kd@beurer.de

youw8
DASHBOARD YOUR LIFE

1



Click on
www.beurer-connect.com

- (TR) Web sitemizi ziyaret edin
(RUS) Посетите нашу страницу в Интернете
(PL) Odwiedź naszą witrynę internetową
(NL) Bezoek onze website

(P) Visite o nosso site na Internet


2



Dilinizi seçin ve kaydolun
Выберите язык системы и зарегистрируйтесь
Wybierz język i zarejestruj się
Selecteer de gewenste taal en registreer uzelf

Selecione o seu idioma e registre-se

3

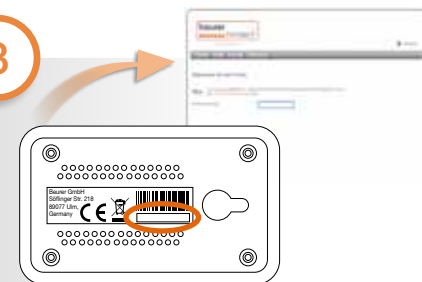

Vericinin etkinleştirme kodunu girin
Введите код активации **радиоустройства**
Wprowadź kod aktywacyjny **modułu radiowego**
Activeringscode van de **draadloze**
ontvanger instellen
Introduzir o código de ativação da **caixa**
de radiotransmissão

4



- (TR) Vericiyi yönlendiriciye ve şebeke kablosuna bağlayın
(RUS) Подключите радиоустройство к маршрутизатору и сетевому кабелю
(PL) Podłącz kabel sieciowy i router do modułu radiowego
(NL) Sluit de draadloze ontvanger aan op de router en het netsnoer
(P) Ligue a box radioelétrica ao router

5



Terazi seçeneğini belirleyin ve etkinleştirme kodunu girin
Выберите **весы** и введите код активации
Wybierz **wagę** i wprowadź kod aktywacyjny
Selecteer de **weegschaal** en voer de activeringscode in
Selecione a **balança** e introduza o código de ativação

6



Vücut değerlerini girin
Введите физические показатели
Wprowadź parametry ciała
Voer de lichaamswaarden in
Introduza os valores corporais

7



- (TR) Pilleri takın, teraziyi zemine yerleştirin
(RUS) Установите батарейки, поставьте весы на пол
(PL) Włóż baterie, ustaw wagę na podłożu
(NL) Plaats de batterijen, plaats de weegschaal op de vloer
(P) Colocar as pilhas, coloque a balança no chão

8



Çıplak ayakla tartılın → Ağırlığı kendi profilinize atayın
При измерении веса снимайте обувь
→ Внесите данные веса в учетную запись
Pomiaru wagi należy dokonywać boso
→ Przypisanie wagi do profilu
Wegen op blote voeten → Koppel het gewicht aan uw profiel
Pese-se descalço → Associe o peso à sua conta

9



Click on App Store/ Play Store

Uygulama ücretsiz üzerinden indirilebilir
Загрузите бесплатно приложение
Darmowe aplikacje do pobrania
De app kan gratis worden gedownload
Aplicação para download gratuito

CE

752.994-1212