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## **SLIMMING BELT**

### **Slimming Belt for Vibratory Massage**

---

**Manual**  
Model B2-12

---

**Thank you for buying vibratory massaging belt  
BODYKRAFT B2-12.**

**Please read the following manual carefully.**

---

Nowadays low motor activity of a man, unbalanced diet, stress, lack of time for physical exercise – all that inevitably leads to functional changes in the organism, such as - irregularities of life activity, premature aging. Massage is scientifically based, time-tested, most physiological for a man's body health-improving way. It can be used extensively by everybody, regardless of age, first of all as a hygienic method of everyday body care, as a part of healthy lifestyle.

Hygienic massage is used to maintain normal functional state of the organism, for prevention of different diseases, and restoration of capacity for work after considerable mental and work load, as a rehabilitation measure of organism restoration after some diseases and traumas.

Massage is especially effective if it's used combined with morning exercises, different sports – jogging, fitness, going to the gym, or combined with bathing procedures. With the help of massage you can get rid of fatigue, restore strength, improve appearance, slow down the aging process, and get back cheerfulness and good spirits.

## PURPOSE AND FUNCTIONS OF A VIBRATORY MASSAGING BELT

### 1.1. PURPOSE OF VIBRATORY MASSAGING BELT

BODYKRAFT vibratory massaging belt is a versatile medical and prophylactic mean of maintaining high standard of organism's functional state, prophylactics of disease, organism's rehabilitation during the process of recovery, improvement of general well-being, appearance and figure. For medical purposes the massager can be used combined with other procedures, recommended by the doctor and as addition to a professional massage.

- Lower back massage
- Calves massage
- Thigh massage
- Waist massage
- Back massage

### 1.2. INDICATIONS FOR USE

*Massager is especially effective for:*

- functional abnormalities of the nervous system, over fatigue, insomnia, heightened irritability and fatigability
- backaches, which occur after significant physical activity and degenerative dystrophic changes of the spine, which are usually a result of sedentary life-style and age-related changes (radiculitis, osteochondrosis)
- muscle atrophy, myositis, joint stiffness, which appear as a result of chronic inflammatory processes, paresis, arthritis and trauma.
- blood flow abnormalities, which origin in venous and capillary congestion and manifest themselves as arms and legs numbness,

swellings and general sense of heaviness in the limbs.

- preparing muscles for physical activity, and stress relieving after it (when doing sports, etc.)
- overweight problem and cellulites, and low skin-muscular tone after pregnancy and childbirth, after quick weight loss.

### **1.3. RESULT OF USE**

- stimulates blood flow, lymph drainage and tissue trophism.
- improves metabolism, vessels flexibility, state of cutaneous covering;
- activates restorative processes in the organism, eliminates stagnant effects;
- enables weight loss and helps to fight with external manifestations of cellulites;
- heightens muscle tone and joints mobility;
- normalizes work of internal organs, boosts immune system;
- helps to return to shape after pregnancy and childbirth;
- prevents flabby skin from appearing after a quick weight loss;
- relieves muscle ache, removes fatigue, heightens work efficiency.


## **TECHNICAL CHARACTERISTICS AND COMPLETE SET OF THE DEVICE**

### **2.1 SPECIAL FEATURES AND ADVANTAGES OF THE DEVICE**

- This model of massaging belt is practical, compact, portable and versatile (as it provides both general hygienic massage of all body parts and figure corrective massage);
- In complete set of the belt there is a stylish sports bag, convenient to carrying and storage. (You can take it with to the gym, to work, cottage, and business trip);
- Compatibility (it's possible to conduct massage without changing your usual schedule, combining it with household chores, gymnastics, active and passive rest.);
- Symmetrical double-point massage, large area of contact of massag-

- ing surface with the body, control of intensity and vibration mode;
- High speed rotating vibration, even distribution of vibrations to tissues, remote control with liquid crystal display and timer
  - Secure fastening of the belt makes regular massage of entire body comfortable, pleasant and easy;
  - Availability of the extender allows people of any constitution use it.

## 2.2 MODEL DESCRIPTION

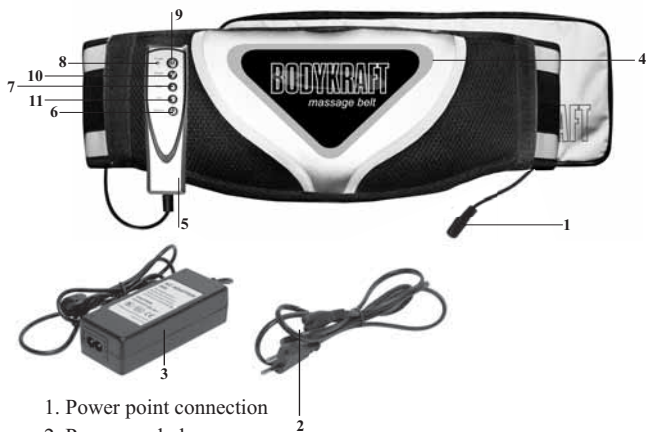1. Power point connection
2. Power cord plug
3. AC Adaptor
4. Vibratory massaging belt
5. Remote control
6. DIRECT button – choice of the direction of massage (clockwise or counterclockwise)
7. HIGH button – massage intensity increase
8. Belt power indicator
9. POWER button – Turning the belt On/Off
10. AUTO/MANUAL button – choice of work mode
11. LOW button – work intensity decrease

## 2.3 TECHNICAL CHARACTERISTICS

Voltage	220-240 V, 50 HZ
Max power	48 W
Number of modes	2
Number of speed levels	5
Automatic shutoff in	10 min
Remote control	+

MANUFACTURER RESERVES THE RIGHT TO INTRODUCE MODIFICATIONS TO THE COMPLETE SET AND TECHNICAL CHARACTERISTICS OF THE DEVICE IN THE PROCESS OF IMPROVEMENT OF IT'S PRODUCT WITHOUT NOTIFICATION ABOUT THOSE MODIFICATIONS.

This product qualifies all required European and Russian standards of hygiene and safety. Made on the order and under supervision of BODYKRAFT LLC, USA. Made in China

## 2.4 COMPLETE SET

- Vibratory massaging belt – 1 p
- Belt extender -1 p
- Power Adaptor -1 p
- Remote Control -1 p
- Bag for transport and storage -1 p

## SAFETY MEASURES

### 3.1 SAFE PRACTICE REQUIREMENTS

- Before use please read the manual carefully.
- Please keep the manual, warranty and, if possible, gift box.
- Before turning on make sure that mains voltage in your house corresponds with voltage on the label of the device.
- Do not keep the device close to heat sources, sunlight and moisture, do not immerse device into the water.
- Do not operate the device with its frame cover open.
- Unplug the device when you're not using it or before cleaning it.
- Do not use the device if you find defects of the cord, plug or other parts.

- Do not try to repair the device yourself.
- Use massaging belt no longer than 20 minutes at a time, and avoid continuous massage of one area for a long period of time.
- Make sure - that children only use the belt under adult supervision.
- Avoid excessive pressure on the device lest you should break it.

### **ATTENTION!**

1. Massage should be pleasant and comfortable, after it you can expect a feeling of pleasant warmth and relaxation in the massaged area. In the beginning of use, during first massage sessions untrained people may feel slight discomfort, which passes in time when massage is regular.
2. If using BODYKRAFT belt causes you pain or strong discomfort, then before using it for the second time consult your doctor.

### **3.2 CONTRA-INDICATIONS**

- Acute feverish states (high temperature).
- Hemorrhaging and predisposition to it.
- Various skin conditions suppurative processes.
- Any acute inflammations of blood and lymphatic vessels, thrombosis, pronounced varicose.
- Mental disorders with excessive excitation states.
- Serious heart, lungs, and other internal organs conditions.
- Acute respiratory diseases.
- Bowels disorder (nausea, vomiting, and liquid stool).

### **WARNINGS**

1. Do not use BODYKRAFT belt to massage traumatized, swollen, inflamed areas or wounds.
2. Do not use the device's high speed to massage your head.
3. Do not massage the area around thyroid gland (just

under Adam's apple).

4. Not recommended for women to use massaging belt while their period and pregnancy.

## OPERATION OF VIBRATORY MASSAGING BELT

### 4.1 HYGIENIC REQUIREMENTS OF MESSAGE CONDUCT

The room where massage is conducted needs to be ventilated.

Massage can be conducted through thin cotton or woolen clean undergarment or cloth, it's necessary to keep the skin clean.

Conducting massage right after ample meal is not recommended, you should wait for 1.5-2 hours.


Massage can be conducted at any time of day and night, in any position that suits you – standing, sitting, lying.

After the massage a brief rest is recommended (10-15 minutes), you can take a shower and apply nourishing or anti-cellulites cream.


### 4.2 WORK ORDER AND MESSAGE CONDUCT SEQUENCE

#### *WORK SEQUENCE*


1. Connect the belt to the power adaptor and plug it in.
2. Put the belt on the necessary part of the body, secure it, and take a comfortable position
3. On remote control (5) push “POWER” button (9), the belt will start working in AUTO mode – indicator (8) will light up green.

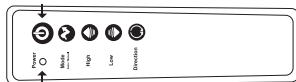
4. “Work speed can be changed with “HIGH” (6) or “LOW” (11) buttons

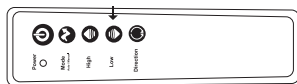5. Work mode can be changed by pushing the button “AUTO/MANUAL” (10)

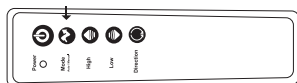6. Direction of massage in a Manual mode can be changed by pushing “DIRECT” button.
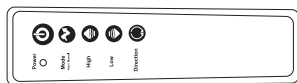
### ***SEQUENCE OF MASSAGE CONDUCT***

Start general hygienic massage session from the back (3-5 minutes) standing or sitting. In doing so – you can secure the belt horizontally or obliquely. Gradually move the belt from upper back and neck to the lower back and pelvis area.

Later conduct massage of the back and side surface of legs – thighs and calves (1-2 minutes), then – the front parts of thighs and stomach. Feet massage can be conducted sitting, without securing the massaging belt, but holding it with your hands, resting feet on it (1-2 minutes)

After that for a minute conduct massage of big pectoral muscles, although avoiding a direct massage of breasts for women.

General body massage can be conducted lying (on the floor, sofa, and bed). You don't have to secure the belt but fold it, using like a massage cushion, resting on it and moving to different body

parts, in the same sequence as the massage standing and sitting.

For a figure corrective massage and fighting the cellulite it's recommended to secure the belt tightly on a local area of the body with large fat deposits (stomach, hips, and buttocks) and conduct massage for 10-15 minutes on max speed on every area. You can finish the massage up with a general body massage.

### ***RECOMMENDATIONS FOR USE***

To conduct a hygienic massage for improving a person's well-being, prevention of premature skin aging, raising the tone and prophylaxis of diseases – a daily conduct of the massage as an addition or substitute of the morning exercises. Back and feet massage is especially effective for people who have to spend a long time in a static position.

Conducting a vibratory massage for a long period of time on a low speed (15-20 minutes) improves blood flow, stimulates muscle relaxation and relieves fatigue, at the same time – a brief high speed massage (4-5 minutes) renders tonic influence, speeds up metabolism, raises muscle tone.

The recommendation is – to start with low speeds and short time, gradually, as body adjusts, speed and time can be increased.

You should measure the force and duration of the massaging belt use, noting your general well-being and subjective sensations, which are characterized by pleasant warmth appearing in the massaged area, after which you should move the belt to the next area.

When doing sports, massage will help to warm up muscles before physical activity, enhance workability, to prevent possibility of sports trauma, relieve fatigue, spasms and remove lactic acid from the muscles after the exercises.

For a figure corrective massage – the best results will be achieved if massage conduct is combined with a healthy diet and exercises (a course no less than 3 weeks with 1 week break). A

norm is – 1-2 massage sessions a day 15 minutes each (morning and evening) max speed, specifically on the problem areas – stomach, waist, and hips. After massage other additional actions can be taken (hot pack, anti-cellulite creams). 3 sessions a day is an allowed maximum, with 2, 5 - 3 hour breaks.

A good result after the massage conduct is – skin reddening, tingling and slight burning sensation, which is a result of the rush of blood to the massaged area and increased metabolism in fat tissue.

### **CLEANING AND ATTENDANCE**

Before cleaning the belt unplug it. If the massaging belt gets dirty, clean it with a cloth, dampened in a soft detergent. Do not use a solvent, benzine and other such agents. Do not use metal brushes and abrasives for cleaning the belt.

### **QUALITY WARRANTY**

This device has a long operation period, on the condition that all the rules of this manual have been followed. Warranty for BODYKRAFT B2-12 massaging belt is 12 months since the day of purchase and covers the defects, which originated in the process of manufacturing. Following recommendations and guidelines, written in the user manual will help you to avoid problems, when operating and servicing the device.

Manufacturer is not responsible for device's defects if service centre proves that malfunctions started after buyer received the device as a result of violation of user manual rules, transportation or storage rules, as a result of third parties' action, force-majeure (fire, natural disaster, etc.); domestic insects or rodents getting inside, exposure to other outside factors and as a result of major violations of technical specifications, which were defined in the user manual, including unstable parameters of the mains, determined by State Standard 13109-97.