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**Thank you for selecting massager BODYKRAFT M-43.
Please, read this manual carefully for correct use. Thank you.**

Low physical activity, unbalanced food, stresses, short time for regular sport leads to functional changes in the organism, and destroy vital functions. Massage is scientifically proved, checked up by long-term practice, the most physiologic for a human body improving means. Independent on age it can be used by each person. First of all, massage will be very helpful as hygienic means in daily care of the body, as a component of a healthy way of life.

Hygienic massage is applied to preservation of a normal functional condition of an organism, preventive maintenance of various diseases, restoration of working capacity after con-

siderable intellectual and physical activities, as a rehabilitation measure on restoration of an organism after some diseases and traumas.

The massage spent in addition to morning gymnastics, playing sports: run, fitness. By means of massage it is possible to take off weariness, to restore forces, to improve appearance, to slow down ageing processes, to return vivacity and good mood.

APPLICATION OF MASSAGER

Massager BODYKRAFT is intended for knead massage of different body parts with four moving adjustable massage heads. Application of the massager is based on dosed mechanical impact of friction and pressure applied directly to the body surface which imitates methods of classical massage such as kneading and rubbing. As a result many various sensitive nerve tergums (receptors) located in the skin, sinews, bands, fascias, muscles, vessels and internal organs are irritated. Impulse flow arising during irritation reaches brain-cortex where all signals are synthesized. This causes body response which is manifested in various positive functional changes in organs and systems such as skin, subcutaneous fat, blood and lymph flow; breathing function and respiratory metabolism; excretion and metabolism and homeostasis; central nervous system and neuromuscular apparatus.

Physical sense of kneading and rubbing is maceration and disintegration of paraplasmic indifferent tissue layers, improving of blood and lymph circulation. Along with it metabolism efficiency increases, contractile function of muscle fibers, their elasticity and mobility improve, muscle lassitude decreases or even eliminates, muscle workability grows. Intensive kneading and rubbing along main nerve trunks and at nerve-endings at the body surface cause decrease of cortex excitability.


FUNCTIONS OF MASSAGER BODYKRAFT

1. Stimulates blood circulation, lymph drainage, tissue trophic;
2. Improves metabolism, vessel flexibility, skin condition;
3. Activates cell functions and cellular secretion, eliminates congestive phenomena;
4. Improves skin and muscles tonus and joints motility;
5. Stabilizes internal organs function, improves immune status;
6. Cures muscle pain, eliminates fatigue cumulation, improves working ability;
7. Promotes relaxation and decreases stress, improved wounds, fractures and trauma healing;
8. Improves condition in chronic fatigue, depression, stress and insomnia;
9. Accelerates restoration after bruises, trauma and fractures.

EFFECT FROM APPLICATION

- Stimulates blood circulation
- Improves a metabolism, flexibility of vessels, an integument condition
- Makes active regenerative processes in an organism, eliminates developments of stagnation
- Raises a muscular tone and mobility of joints
- Normalises work of an internal, strengthens immune system
- Allows to struggle effectively with external displays of a cellulitis
- Prevents flabby skin after fast growing thin
- Facilitates pains in muscles, removes feeling of weariness, raises working capacity
- Accelerates regeneration of the damaged fabrics that promotes restoration after bruises, traumas and crises
- Improves a condition at chronic weariness, depression, stress and a sleeplessness
- Reduces pains at arthritises, pojasnichno-kresttsovom a radiculitis.

TECHNICAL DATA AND COMPLETE SET

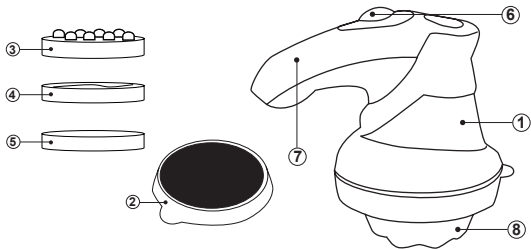1. Massager
2. Protective cap
3. Nozzle with rollers for deep stimulation
4. Nozzle with a wavy surface for strengthening of effect of massage.
5. Nozzle with a smooth surface for soft influence.
6. Regulator of operating modes
7. The masseur handle
8. Massage head.

FEATURES AND ADVANTAGES OF THE PRODUCT

- Compactness and ease.
- The ergonomic design вибромассажера BODYKRAFT provides a maximum of simplicity and convenience of application.
- The convenient form прорезиненной handles prevents device sliding in a hand and allows to make massage of all body at any angle.
- The device case is executed from the non-polluting qualitative practical plastic applied in medicine.
- Rotating massage head, high speed of rotation.
- High intensity of vibration, uniform transfer of fluctuations to fabrics.
- Smooth adjustment provides a wide range of a choice of modes of vibration and intensity of influence.
- Demountable massage nozzles for carrying out of various kinds of massage.
- Protective cap for safe use.

TECHNICAL DATA

Type	Massager
Model №	M-43
Voltage	220-240 V., 50 Hz
Power	25 W
Adjustment of intensity of massage	+
Attachments	3
Protective cap	+
Smooth adjustment of intensity	+
The general time of continuous use	20 min.

THE MANUFACTURER HAS THE RIGHT TO INTRODUCE MODIFICATIONS TO DESIGN, COMPLETE SET, AS WELL AS TECHNICAL DATA OF THE DEVICE IN THE COURSE OF CONTINUOUS IMPROVEMENT OF THE PRODUCTS WITHOUT ADDITIONAL NOTIFYING OF SUCH MODIFICATIONS.

THIS PRODUCTS CONFORMS WITH ALL EUROPEAN AND RUSSIAN SAFETY AND HYGIENIC STANDARDS APPLICABLE.

COMPLETE SET

Manual – 1 pc.

Massager – 1 pc.

Cord – 1 pc.

Massage attachment – 3 pcs.

Protective cap – 2 pcs.

SAFETY USE

- Read this manual carefully before using masager

BODYKRAFT

- Please, keep this manual, warranty certificate and package, if possible.
- The instrument is designed for home use only.
- Keep the instrument protected from heating sources, high temperature, sunlight and moisture.
- Don't immerse the instrument into water.
- Don't use this instrument in basins, bathrooms and other places with high moisture.
- Don't use this instrument with open body cover.
- Supervise children using masager BODYKRAFT any time.

ATTENTION!

- 1. MASSAGE SHOULD GIVE COMFORT AND PLEASANT FEELING. IF USING MASSAGER BODYKRAFT YOU FEEL PAIN OR DISCOMFORT CONSULT YOUR PHYSICIAN BEFORE NEXT SESSION.**
- 2. DON'T USE MASSAGER ON DAMAGED OR IRRITATED SKIN.**

CONTRA-INDICATIONS TO CARRYING OUT MASSAGE

- Sharp inflammatory processes and feverish conditions (heat).
- Bleedings and predisposition to them.
- Tumours of various localisation
- Some diseases of a skin, purulent processes.
- Inflammations of lymph nodes and vessels
- Thrombosis or the acute inflammation of veins expressed varix dilatation.
- Mental diseases with excessive excitation.
- Sharp respiratory diseases.
- Infringements of work of intestines (a nausea, vomiting, a liquid chair).

WARNING

- Application of Massager BODYKRAFT is prohibited for persons with implanted stimulators.
- Don't use massager BODYKRAFT for massage of injured, swollen or irritated areas or wounds without your physician's advise.

- Massage with application of this instrument is prohibited for children younger 12 years.
- Application of the instrument is not recommended for pregnant women.

APPLICATION AND PROCEDURE OF MASSAGER

1. Turn the cover, then put the 3 pcs AAA batteries into the battery's box according to the symbol «+» and «-» on the unit.
2. Press on the top of cover, the massager starts working.

HYGIENE OF MASSAGE

Massage room should be ventilated before the session starts. Body and massage surfaces should be clean. You may perform massage through thin cotton or wool underwear. It is

not recommended to perform massage session immediately after the meal, you should wait about 2 hours.

You may perform massage at any time during the day, in any comfortable posture – standing, lying, sitting. Make effort to relax muscles in the area of massage.

Short passive relax (10 – 15 min) is recommended after massage session, you may take a shower and apply moistening and nourishing cream on your skin.

RECOMMENDATIONS FOR APPLICATION

For hygienic massage to improve well-being, to protect skin and body from ageing, to improve tonus and disease prevention. Massage is recommended as a daily additional technique or substitution to morning exercises.

Recovery massage helps to recover lassitude after sport exercises, improve workability, prepare yourself for physical stress, warm up muscles and prevent sport trauma. Usually 1 – 2 sessions are performed daily, it is allowed to make 3 massage sessions during one day with 2.5 – 3 h intervals. Daily sessions are most effective, intervals longer than 2 days are not desirable.

Massage should not cause increase of pain intensity. You may expect feeling of heat, comfort, relaxation in the massaged area, improvement of general well-being, of joints flexibility, of

respiratory easiness and freedom.

MASSAGE PROCEDURE

Start session of general hygienic massage from shoulder (5 – 10 minutes) in sitting, laying or standing posture. Mount massager on support and press its massage surface to the necessary body area. Any horizontal or vertical surfaces may be used as support, including walls, cupboards, armchair, chair, coach back; bed or floor. Move the instrument gradually, massage body sides (each place 1 – 2 minutes). Further massage rear leg surfaces – thighs and calves (1 – 2 minutes), and then – front and side parts of thighs and belly (1 – 2 minutes).

Feet massage may be performed in sitting posture with the pil-low laying on the floor or in the laying posture with legs bended (1 – 2 minutes). After that massage greater pectoral muscles during 1 minute avoiding direct strong massage of mammary glands.

Reverse function should be used for periodical alteration of rotation direction of massage rolls. Alteration is used according to the following principle: half session – in one direction, the other half – in the opposite direction.

It is recommended to start massage from mild massage and with subsequent gradual increase of depth and strength and also duration of the session.

CLEANING AND MAINTENANCE

If contaminated, clean it with tissue moistened with mild cleaning agent. Don't use solvents, benzene and similar agents. Don't use metal brushes or abrasives for cleaning.

WARRANTY

This instrument is intended for long-term use if requirements of this manual are observed. Warranty on BODYKRAFT massager is 12 months from the date of purchase and covers defects arising from manufacturing. Following recommendations and instructions of this manual (customer's manual) ensures absence of trouble during use and maintenance of this instrument.

Please, note that this device is intended for home use only, not for business purposes. Application of this instrument for other purposes not mentioned above is violation of rules of correct use of the device.

Manufacturer is not liable for defects if the service department proved that those defects arose after the instrument was transferred into customer's possession as a result of violation of rules for use, transportation, storage, third party's actions, circumstances beyond control (fire, natural disaster, etc.); home insects or gnawers impact, influence of foreign factors as well as

significant violation of technical requirements mentioned in the manual, including unstable power supply, established in GOST 3109-97.