OMRON

- (EN) Instruction Manual
- (IT) Istruzioni per l'uso
- (DE) Gebrauchsanweisung
- (FR) Mode d'emploi
- (NL) Gebruiksaanwijzing
- (ES) Manual de instrucciones
- (RU) Руководство по эксплуатации
- (AR) كتيب الإرشادات
- (BG) Наръчник с инструкции
- (CS) Návod k obsluze
- (DA) Brugsanvisning
- (EL) Εγχειρίδιο Οδηγιών
- (ET) Kasutusjuhend
- (FI) Käyttöohje
- (HR) Upute za korištenje
- (HU) Használati útmutató
- (LT) Naudojimo instrukcija
- (LV) Lietošanas instrukcija
- (NO) Bruksanvisning
- (PL) Instrukcja obsługi
- (PT) Manual de Instruções
- (RO) Manual de utilizare
- (SL) Navodila za uporabo
- (SK) Návod na obsluhu
- (SV) Instruktionsmanual
- (TR) Kullanma Kılavuzu

HN286 Digital Personal Scale


A Good Sense of Health

IM-HN-286-E-01-05/2010
5329447-7A

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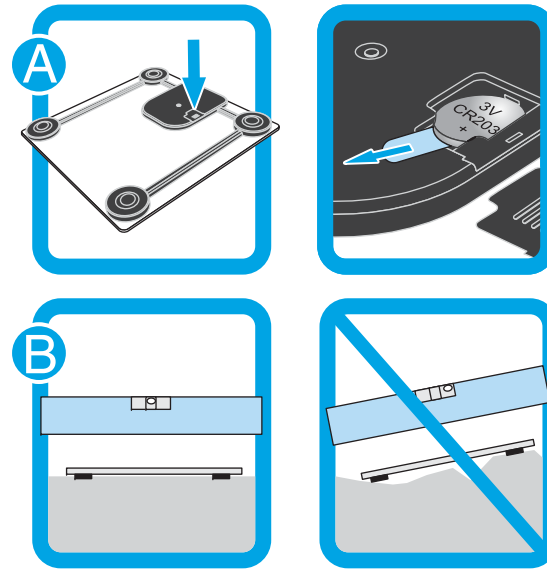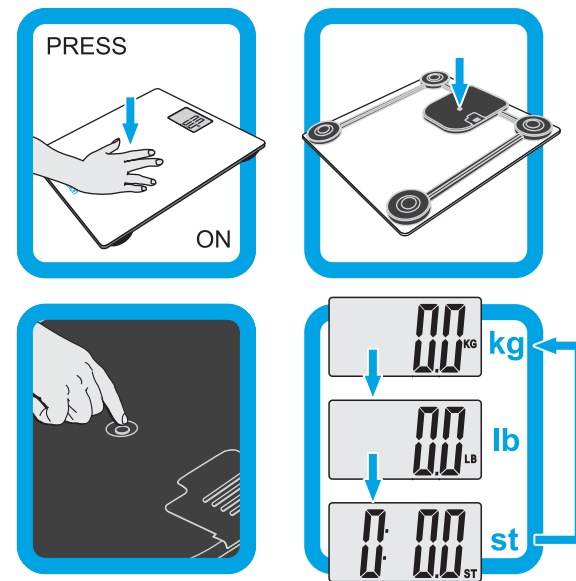
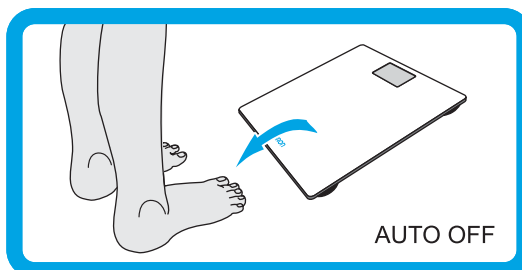
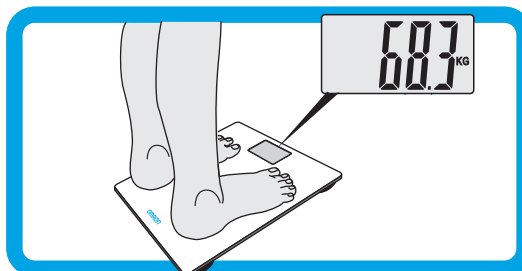
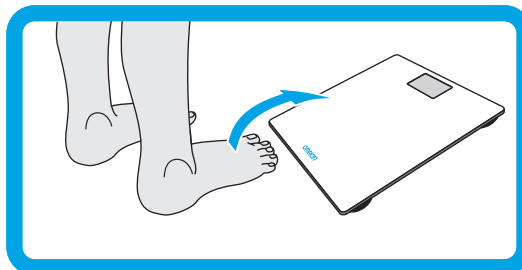
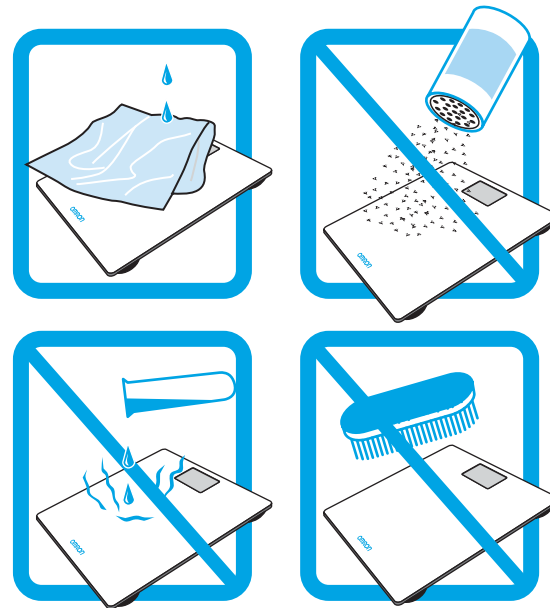
Symbols and definitions are as follows:

Warning: Improper use may result in possible death or serious injury.

Caution: Improper use may result in injury or property damage.

Warning:

- Never start weight reduction or exercise therapy without the instructions of a doctor or specialist.
- Do not use the scale on slippery surfaces, such as a wet floor.
- Keep the scale out of the reach of young children.
- Do not jump onto the scale, or hop up and down on the scale.
- Do not use this scale when your body and/or feet are wet, such as after taking a bath.
- Do not step on the edge or display area of the scale.
- People with disabilities, or who are physically frail, should always be assisted by another person when using this scale.
- Use a handrail or so when stepping on the scale.
- If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Consult a doctor immediately.
- Do not put your fingers into any slot or inside the scale.
- Use indoors.
- Do not use for purposes other than taking weight.

Caution:

- Approved for human use only.
- The scale is slippery glass, be carefully of falling.
- Always clean the scale before use when using the unit with people with a skin disease on foot.


- Do not store in locations that exposed to chemicals or corrosive vapours.
- Please read this instruction manual carefully before use and for further information on the individual functions.
- Do not place the scale against the wall. May cause injured and malfunction.
- Do not disassemble, repair, or remodel the scale.
- Do not use battery not specified for this scale. Do not insert the battery with the polarities in the wrong direction.
- Replace worn battery with a new one immediately.
- Do not dispose of batteries in fire.
- If battery fluid should get on your skin or clothing, immediately rinse with plenty of clean water.
- Remove the battery from this scale when you are not going to use it for a long period of time (approximately three months or more).
- This scale is intended for home use only. It is not intended for use in hospitals or other medical facilities. It does not support the quality standards required for professional use.

General Advice:

- For correct measurements:
 - Please store this scale on level ground.
 - Do not place this scale on a cushioned floor surface such as on a carpet or a mat.
 - Do not store any objects on the platform.

EN Messages
 IT Messaggi
 DE Meldungen
 FR Messages
 NL Meldingen
 ES Mensajes
 RU Индикация на дисплее
 AR الرسائل
 BG Съобщения
 CS Hlášení
 DA Meddelelser
 EL Μηνύματα
 ET Teated
 FI Huomautuksia
 HR Poruke
 HU Üzenetek
 LT Užrašai ekrane
 LV Ziņojumi
 NO Meldinger
 PL Komunikaty
 PT Mensagens
 RO Mesaje
 SL Sporočila
 SK Indikácie na displeji
 SV Meddelanden
 TR Mesajlar

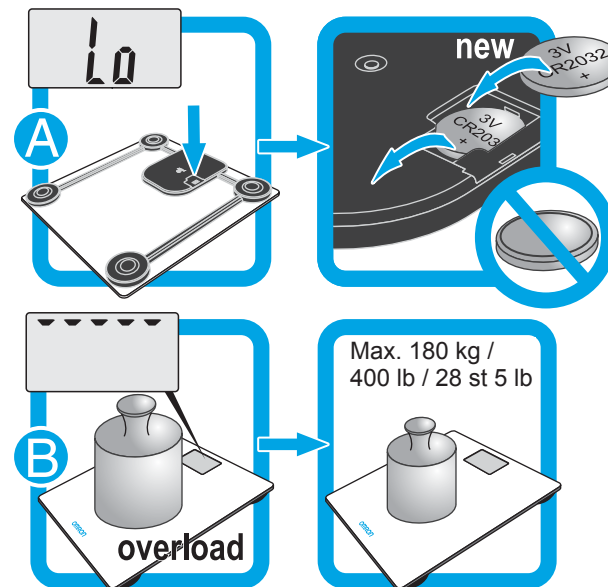Name	HN286 (HN-286-E)
Product	Digital Personal Scale
Display	0 ~ 180 kg 0.1 kg 0~400 lb: 0.2 lb 0~28 st 5 lb: 0.2 lb
Weight Accuracy	5.0 to 60.0 kg: ± 0.4 kg 60.1 to 100.0 kg: ± 0.6 kg 100.1 to 180.0 kg: ± 1%
Power Supply	One CR2032 battery (test battery is included)
Battery Life	Approximately 1 year (when scale is used for 4 measurements a day)
Operating Temperature / Humidity	+5°C to +35°C, 30% to 85% RH
Storage Temperature/ Humidity/Air Pressure	-20°C to +60°C, 10% to 95% RH, 700 hPa - 1060 hPa
Weight	Approximately 1.5 kg (battery included)
External Dimensions	Approximately 300 (W) × 18 (H) × 300 (D) mm
Package Contents	Digital Personal Scale, one CR2032 battery for test use of the scale, instruction manual

Note: Subject to technical modification without prior notice.

This device fulfills the provisions of the EMC Directive 2004/108/EC.

OMRON HEALTHCARE EUROPE B.V. guarantees this product for 2 years after date of purchase. The guarantee does not cover battery, packaging and/or damages of any kind due to misuse (such as dropping or physical misuse) caused by the user. Claimed products will only be replaced when returned together with the original invoice / cash ticket.

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OMRON HEALTHCARE Co., Ltd.
OMRON HEALTHCARE UK LTD.
OMRON Santé France SAS