

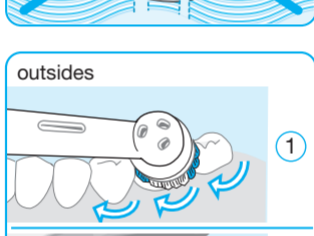
# Oral-B®


VITALITY powered by BRAUN

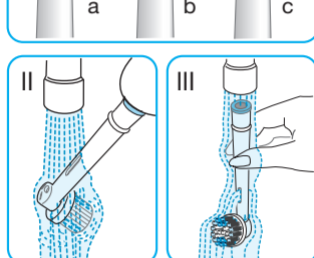
Type 3757

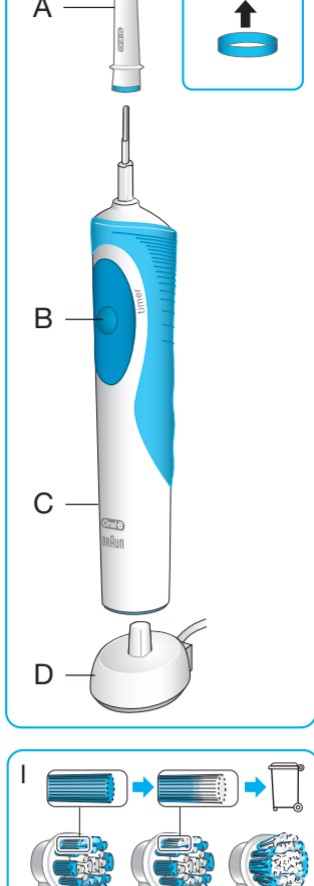
outsides

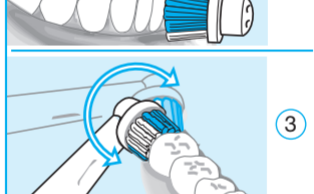insides

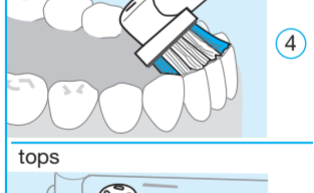tops

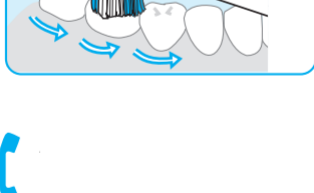- Internet:**  
[www.oralb.com](http://www.oralb.com)  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)
- Braun GmbH  
 Frankfurter Straße 145  
 61476 Kronberg/Germany
- D12.013  
 D12.513  
 D12.513W  
 D12.513P  
 D12.513S  
 D12.523

## Deutsch

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### English

Your Oral-B Vitality toothbrush has been developed with dental professionals to effectively remove plaque from all areas of your teeth.

Please read the use instructions carefully before first use.

### IMPORTANT

- Periodically check the cord for damage. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre. A damaged or non-functioning unit must no longer be used.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety. In general we recommend to keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance/charger.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

### Description

A Brush head (depending on model/ interchangeable)  
B On/off switch  
C Handle  
D Charging unit

### Specifications

For voltage specifications please refer to the bottom of the charging unit.

### Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet. Place the handle (C) on the charging unit.

- A full charge takes 16 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).

- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible. However for environmental reasons,

Oral-B recommends to unplug the charging unit until the next recharge is required.

### **Battery maintenance**

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

### **Using the toothbrush**

The toothbrush can be used with any leading brand of toothpaste.

- Wet brush head and apply toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (B).
- Slowly guide the brush head from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces (1–5). Brush at least 2 minutes to thoroughly remove plaque. Do not press too hard or scrub.
- Turn off the handle by pressing the on/off switch.


### **Timer (depending on model)**

An incorporated timer memorises the elapsed brushing time. A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

### **Individual brush heads and accessories**

Oral-B offers you a variety of brush heads which are available in stores that sell Oral-B products. Following brush heads are interchangeable and can be used on your Oral-B Vitality handle.

#### Oral-B FlossAction brush head

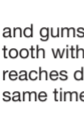The Oral-B FlossAction brush head is recommended for everyday deep cleaning. It features micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque.

#### Oral-B Precision Clean brush head


The Oral-B Precision Clean brush head provides effective and gentle cleaning action to help protect both teeth and gums. The brush head cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

#### Oral-B Sensitive brush head


The Oral-B Sensitive brush head provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums.

#### Oral-B Dual Clean brush head

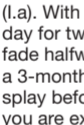The Oral-B Dual Clean brush head offers two types of cleaning motion in one. It combines the oscillating–rotating motion with an additional side-to-side motion. The oscillation motion stimulates gums while removing plaque bacteria. Multi-angled bristles clean between teeth and along the gumline. You can use this brush head either by moving from tooth to tooth or apply the same brushing technique as you do with any manual toothbrush.

#### Oral-B 3DWhite brush head

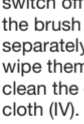For occasional or everyday polishing, please use the Oral-B 3DWhite brush head. It features a dentist-inspired polishing cup to remove surface stains from coffee, tea and tobacco for effective brightening. We do not recommend using the Oral-B 3DWhite brush head with braces.

Dentists recommend replacing any toothbrush every 3 months. Therefore the brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need (I.a). With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period (I.b). If the bristles splay before the colour recedes (I.c), you are exerting too much pressure on teeth and gums.


During the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist.

### **Cleaning**

After use, rinse the brush head for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brush head. Rinse both parts separately under running water (III) and wipe them dry. From time to time, also clean the charging unit with a damp cloth (IV).

Subject to change without notice.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.
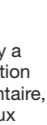
### **Guarantee**

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

### **For UK only:**

This guarantee in no way affects your rights under statutory law.

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