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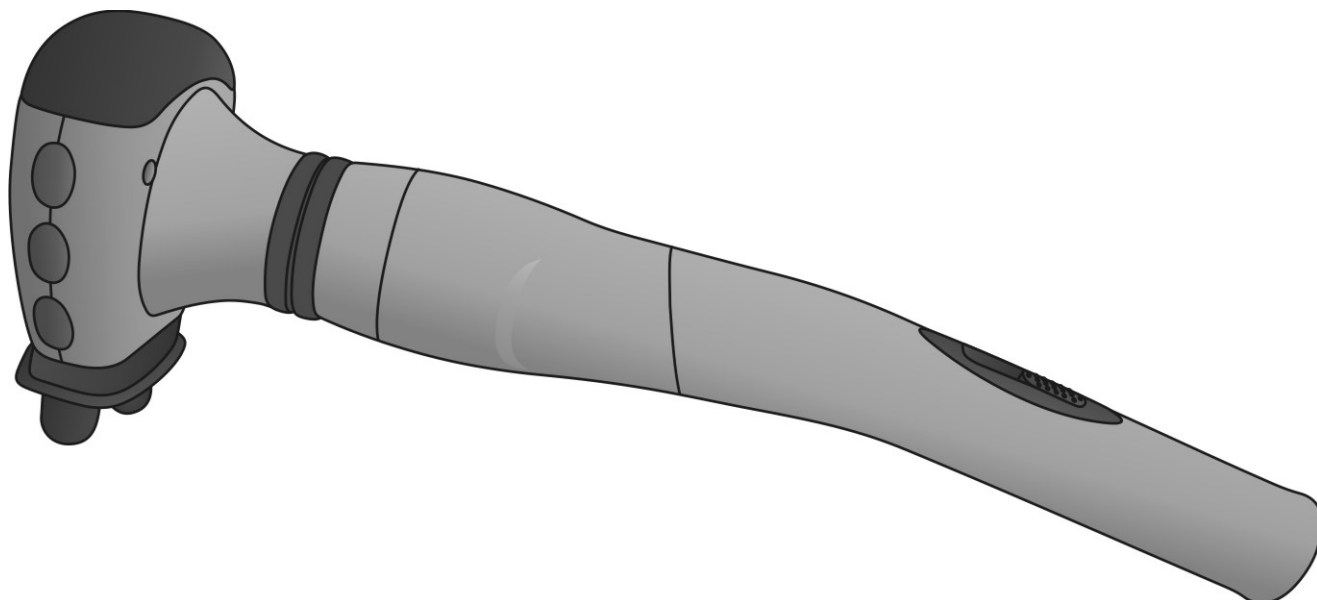
## **GB** INSTRUCTION MANUAL

### **IMPORTANT SAFEGUARDS**

- Please read this instruction manual carefully before use and keep in a safe place for future reference.
- Before the first connecting of the appliance check that voltage indicated on the rating label corresponds to the mains voltage in your home.
- For home use only. Do not use for industrial purposes. Do not use the appliance for any other purposes than described in this instruction manual.
- Do not use outdoors.
- Always unplug the appliance from the power supply before cleaning and when not in use.
- To prevent risk of electric shock and fire, do not immerse the appliance in water or any other liquids. If it has happened DO NOT TOUCH the appliance, unplug it immediately and check in a service center.
- Do not take this product to a bathroom or use near water.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not leave the appliance switched on when not in use.
- Do not use other attachments than those supplied.
- Do not operate after malfunction or cord damage.
- Do not attempt to repair, adjust or replace parts in the appliance. Check and repair the malfunctioning appliance in the nearest service center only.
- Keep the cord away from sharp edges and hot surfaces.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not pull, twist, or wrap the power cord around the appliance.
- Never drop or insert any object into openings on the appliance.
- Do not put operating appliance on any surface, it may cause injury or damage.
- The massager is designed for intermittent use: it can be applied without interruption not longer than **20** minutes, then the motor must be allowed to cool.
- Do not use the massager if you are pregnant, or suffer from thrombosis, open wounds, broken skin and varicose veins. Consult your doctor first.
- Immediately interrupt procedure if you feel discomfort, or if any pain or irritation results. First, obligatory consult your physician.
- The appliance has a hot surface. Insensitive to heat persons should use the appliance carefully.
- If the product has been exposed to temperatures below 0°C for some time it should be kept at room temperature for at least 2 hours before turning it on.
- The manufacturer reserves the right to introduce minor changes into the product design without prior notice, unless such changes influence significantly the product safety, performance, and functions.

### **GENERAL INSTRUCTIONS**

- Use of this massager has a relaxing, soothing and stimulating effect on your muscles. If used regularly for several minutes a day, it can reduce stress, fatigue and muscular pain, and it has anti-cellulite function.
- When using skin or face remedies, this appliance helps to spread the product more evenly and thus enhance its effective skin penetration.
- This appliance is equipped with an infrared lamp. The infrared radiation penetrates into the deep layers of skin causing the pores and vessels to expand and simultaneously enhancing the metabolism process and microcirculation. Using the infrared massager brings relief from muscular pain and regenerates the skin.

### **ACCESSORIES**

#### **ROLLING MESSAGE ATTACHMENT**

- The rolling attachment allows to move the massager on body freely, thus to carry out consecutive transition from one site of body to another, and special ledges on rollers make massage more effective.

#### **ACCUPOINT MESSAGE ATTACHMENT**

- The accupoint attachment is intended for more gentle massage of big sites of body such as back, hips, and legs.

#### **MICRO-FINGERS MESSAGE ATTACHMENT**

- Renders more intensive massage, and daily application of this attachment will allow you to support body muscles in tone and be always in good form. Massage by circular movements, slightly pressing on skin to increase efficiency of procedure.

#### **CONTOUR BODY ATTACHMENT**

- Massage with contour body attachment improves blood circulation and reduces accumulation of salts as well as cellulite efficiently. Use on body directly without clothes.

#### **MULTI-ACCUNODE MESSAGE ATTACHMENT**

- This attachment provides deep and strong penetrating massage to large muscle area for relaxation and soothing muscle pain.

#### **NOTE:**

- To ensure optimum effectiveness during use of any function, move the unit slowly, in either circular or back and forth motion.
- Always begin at a spot that is most distant from the heart. Then gradually work your way closer to the heart, using the motions described.
- Never treat the same spot longer than 3 minutes at a time.

### **ADJUSTMENT**

#### **ROTATION OF MESSAGE HEAD**

- The message head can be rotated downward 180°.

#### **TWISTED HANDLE**

- You can turn the handle on 180° for simplification of access to any sites of a body

### **INSTRUCTIONS FOR USE**

- Plug the appliance to the power supply and slide the function switch to the desired position:

0 – Off.

1 – VIBRATION MESSAGE

- Enclosed attachments allow providing the right massage mode for each part of your body.

- **NOTE:** Although the massager may be applied through clothing, it is most effective when acts directly on the body.

## 2 – INFRA-RED HEATING

- Infrared radiation penetrates deep into muscles, enhancing blood circulation and reducing tension in muscle areas. Use the infrared heat function can also accelerate the restorative processes in the case of colds, inflammation and minor injuries.

## 3 – VIBRATION MASSAGE WITH INFRA-RED HEATING

- Combination of infrared heat and vibration massage intensifies effect of both functions.

### CLEAN AND CARE

- Switch the appliance off and ensure that the massager has completely cooled.
- Wipe over the housing with a soft damp cloth and then with a dry cloth.
- Do not use hard abrasives or cleaners.
- Attachments may be cleaned in soapy water

### STORAGE

- Be sure that the appliance is unplugged and cooled down.
- Complete instructions of section CLEAN AND CARE
- Keep the appliance in a cool, dry place.

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