sinbo

SMR 4229 MASAJ ALETİ

KULLANMA KILAVUZU

TR

EN

DE

FR

ES

NL

RU

HR

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Children being supervised not to play with the appliance.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

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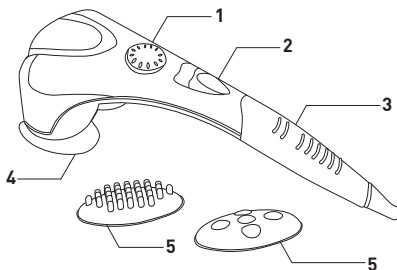
SAFETY NOTICE MATTERS BEFORE USING

- If the patient wants to use this massage machine, please ask for the doctor whether it suits or not.
- Please do not leave when the machine is running especially with the presence of children.
- Please do not cover the machine with any goods when the machine is running.
- Please do not run the machine more than 20 mins.
- Please do not run the machine in the long time, owing to the surface of the machine will be hotter and hotter, after please wait it for cooling to use once more.
- Please do not put this machine on the dry-skin surface.
- This machine only supply the massage non-medical thing, please do not use in the therapy.
- Please do not use this machine before sleeping; it will affect your sleeping.
- Please do not run this machine in the bed.
- Please do not run this machine on the wound.
- The children do not running the machine except to have adult there.
- The surface of the machine will be hotter and hotter when it is running. Please carefully use it under the vibration massage hot compress condition.

USAGE

- Please make sure the pattern of massage stick point at "OFF" before power on.
- Power supply: AC230V
- Adjust the pattern of the massage stick point at "Massaging" or "Massaging and hot compress function"

OPERATION MODES

- Adjust the speed button of the massage stick.

Operation Mode Off

Massaging

Massaging and compress function

- If you running the machine according to the instruction, will make you relax, do not worried about harm to your muscle.
- Only need to move back and forth, each place do not surpass 3 mins. And adjust the place and time of the massage according to your requirement in it order to avoid damage.

- Attaches two kinds of detachable cushion(finger-pressure type and carding brushing type) install it according to the chart.


- Please depend on your enjoy adjust the speed of the vibration.

CHART OF BODY POINT

There are around 650 points in our body from neck to waist. And every point connected with some organs or body condition. Massaging can promote blood circulation, promote metabolism and build up the ability againts disease.

Massager will make you have massive massaging joys and greatly pleasure by easy ways.

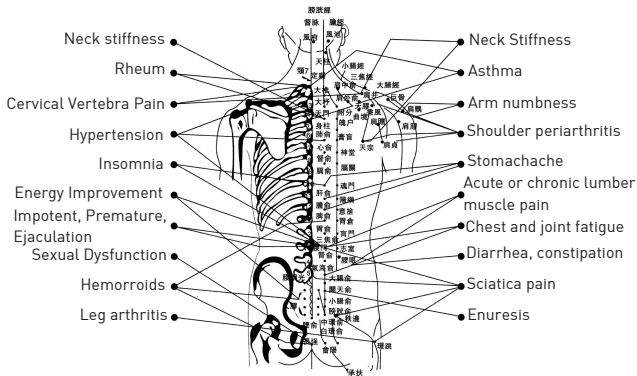THE IMPORTANT NOTICE MATTER

When use this product, should comply for the instruction

Firstly, please read the instruction before running the machine to know danger in the running.

- Please cut off the power after running.

THE IMPORTANT NOTICE MATTER

- Please do not keep this machine in the water or easy to touch water.
- Please do not running this machine in the bathroom.
- This machine only use in the family, not in the business.


Warning: In order to avoid using this machine to create the injury like fire, the electric wire sparking and so on. Please notice following:

- Please do not add other parts in privately when uses this machine.
- This machine does not apply for the children handicapped and patient.
- Please keep this machine away from the fire.
- Please do not running this machine in the footcloth or bedclothes all the time in order to avoid to take fire.
- Please do not use this machine in the sleeping.
- Please do not open the snugly systems when install the detachable cushion.
- Please do not keep close to eyes and sensitive body parts.
- Please do not remove the machine in privately in order to avoid damage the machine.
- Please do not running this machine in the soft cloth.
- Please do not running this machine in the open air.
- Please make sure the power is turned off, after you running.

THIS MACHINE CHARACTERISTIC

- Vibration massage works together with hot pack leisure device, enjoys double from one massage.
- The streamline long handle design it is very easy to take and operate.
- Non-step speed control you can adjust the machine freely.
- Attaches two sets detachable massage parts.
- Intellectualized (Intelligent) time control design safe and comfortable.

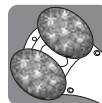CLEANOUT AND MAINTENANCE METHOD

Cleanout

- When the machine body keep the cold and make sure the power will be cut so that you could use dry cloth to clean it.
- Please do not drop the water or liquor into the massage's handle.

Collection

- When you put electrical source out from the socket please handhold the pin by your hands that remember not to drag the electrical wire and keep the machine in the dry place. Please do not enlace the wire together in order to avoid damaging the wire.

JUDGEMENT OF FAULT TROUBLE SHOOTING

Fault	Reason	Fault Remedy
Motor not working	Plug not fixed up Switch is out of position	Fix up the plug well Checking the switch
Working too long	Thermal Protection	Rest for 5 minutes

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United Favour Development Limited
Botek Testing Laboratories Co., Ltd.