**INFRARED MASSAGER**

- Part names**
1. Massage head
  2. Infrared massage surface
  3. 270°-rotation ring
  4. Upper handle
  5. Lower handle
  6. Massage switch
  7. Infrared heat switch
  8. Cord
  9. 90°-bending knob

- Safety instructions**
- Before using the appliance for the first time read this manual carefully.
  - Please keep this manual, the warranty and the packing if possible.
  - The appliance is designed for household use only.
  - Keep the appliance away from heat sources, direct sunlight and humidity.
  - Never immerse the appliance in water.
  - Unplug the appliance with dry hands only.
  - Always unplug the appliance when not in use and before cleaning.
  - Never pull on the cord to unplug the appliance, grip the plug instead.
  - Never allow children to use the appliance unattended.
  - Do not operate the appliance if the cord, the plug or any other part is damaged.
  - Never try to repair the appliance yourself.
  - Do not operate the appliance for more than 20 minutes at a time and avoid treating the same area for a long time.
  - Before turning on the massager make sure that local AC power corresponds to that indicated on the specification plate of the appliance.
  - Do not apply excessive pressure as this may cause damage to the appliance.

- Do not operate the appliance for more than 20 minutes at a time and avoid treating the same area for a long time.
- Before turning on the massager make sure that local AC power corresponds to that indicated on the specification plate of the appliance.
- Do not apply excessive pressure as this may cause damage to the appliance.

- Do not operate the appliance for more than 20 minutes at a time and avoid treating the same area for a long time.
- Before turning on the massager make sure that local AC power corresponds to that indicated on the specification plate of the appliance.
- Do not apply excessive pressure as this may cause damage to the appliance.

- Warning!**
1. If you feel discomfort after using the massager you should consult your doctor before using the appliance again.
  2. Do not use the massager if you have any skin damages or dermatitis.

**Functions of your massager**

- soothes tiredness;
- improves blood circulation;
- activates cell metabolism;
- improves functioning of the internal organs;
- helps to do away with rheumatic and muscular ache;
- helps to relax and relieve muscular tension;
- improves skin condition;
- helps to do away with vertebral ache.

**Infrared effects**

Infrared heat ensures deep stimulation of skin and blood vessels and improves blood circulation. Besides, infrared heat opens the skin's pores and stimulates cell rebuilding.

Infrared heat is applied both for medical and for cosmetic purposes.

**Medical application**

Deep heating helps to remove rheumatic ache and other similar conditions. If you have any of such diseases you should consult your doctor before using the appliance for the first time.

**Cosmetic application**

Use cosmetic creams and lotions to achieve a more profound infrared effect. By intensifying blood circulation, infrared heat improves the state of your skin and makes it healthy and clean.

**Using the massager**

- Turning on**
- This massager is designed for household use only.
  - Connect the cord to an AC outlet.
  - Set the 3-position switch to one of the following modes:
    - 0 - off;
    - 1 - relaxation massage;
    - 2 - intensive massage.
  - To use the infrared massage surface set the infrared heat switch to the utmost position.

**Caution!**

- Connect the plug to a 230 V power source only.
- Do not operate the appliance if the cord or the plug is damaged.
- If the massager does not work properly in any way, unplug it immediately and refer to a qualified technician for repair.

**Vibration massage with infrared heat**  
This appliance allows combining massage with heating. Set the function switch to the "3" position. Vibration massage combined

with infrared heat is most effective when used without any additional attachments. Massage combined with deep heating facilitates muscle relaxation and removes rheumatic ache and other similar conditions.

**Infrared heat**  
Test the temperature of the infrared plate before use. You may also adjust the temperature according to your preference. To reduce the temperature you should periodically set the function switch to the "0" position for a few seconds, then again to the "3" position.

**Vibration massage**  
Six interchangeable attachments are provided for massaging every part of the body: arms, legs, back, shoulders, thighs, neck and face. Place the desired attachment on the massage head. The appliance works most efficiently if you use it without clothing.

1. Acupuncture massage attachment, which provides deep massage to muscles and soothes tiredness.
2. Double rolling massage attachment for scalp massage and improvement of blood circulation.
3. Gentle massage attachment, which provides massage of arms, wrists, legs and shoulders.
4. Shiatsu massage attachment, which soothes muscular tiredness and improves blood circulation.
5. Thump massage attachment, which provides deep massage of the area being treated.
6. Multi-nodes massage attachment, which provides a deep penetrating massage to large muscular areas.

For the sake of convenience, you may bend the massage handle up to 90°. To do so press and hold the bending knob. Release the knob to fix the desired bending angle. This massager features a 270°-rotation ring that allows massaging of hard-to-reach places.

**After use**

It is not recommended to operate the appliance for more than 20 minutes at a time. Wait for at least 20 minutes before using the massager again.

- Turn off the appliance as follows:
1. Set the 3-position massage switch to the "0" position - off.
  2. Set the infrared heat switch to the uppermost position.
  3. Unplug the appliance. Do not pull on the cord; grip the plug instead.

**Cleaning**

- Unplug the massager before cleaning.
- Wipe the outer surface of the massager with a damp cloth. Use mild detergents to remove contamination. Never use metal brushes and abrasives.

**Specifications**





Power requirement: 230 V – 50 Hz  
Maximum power: 13 W  
Recommended working time: 20 minutes

*The manufacturer reserves the right to make changes to the technical characteristics of this device with out prior notification of the consumer.*

**Service life – 3 years**

**Guarantee**  
Details regarding guarantee conditions can be obtained from the dealer from whom the appliance was purchased. The bill of sale or receipt must be produced when making any claim under the terms of this guarantee.

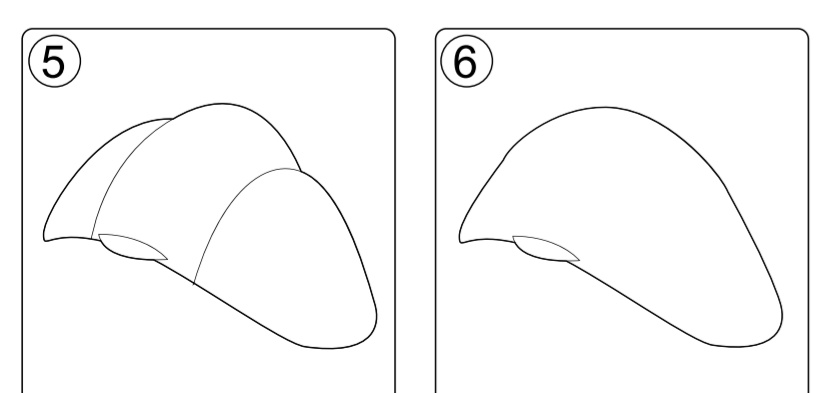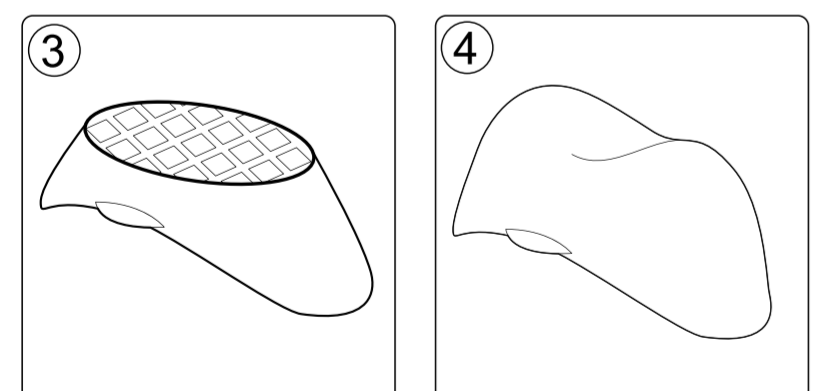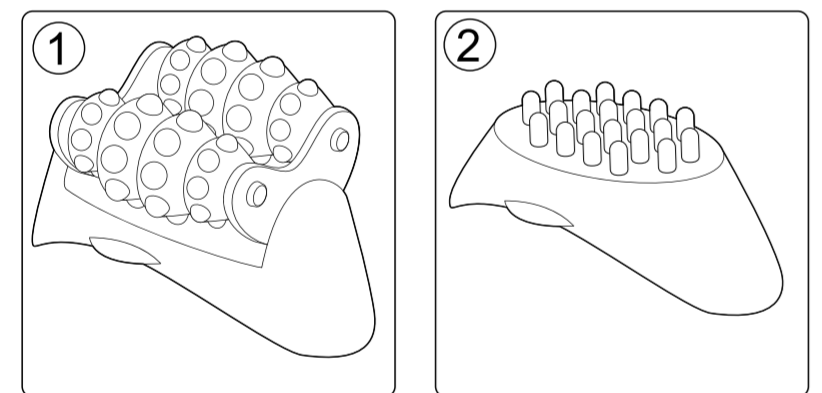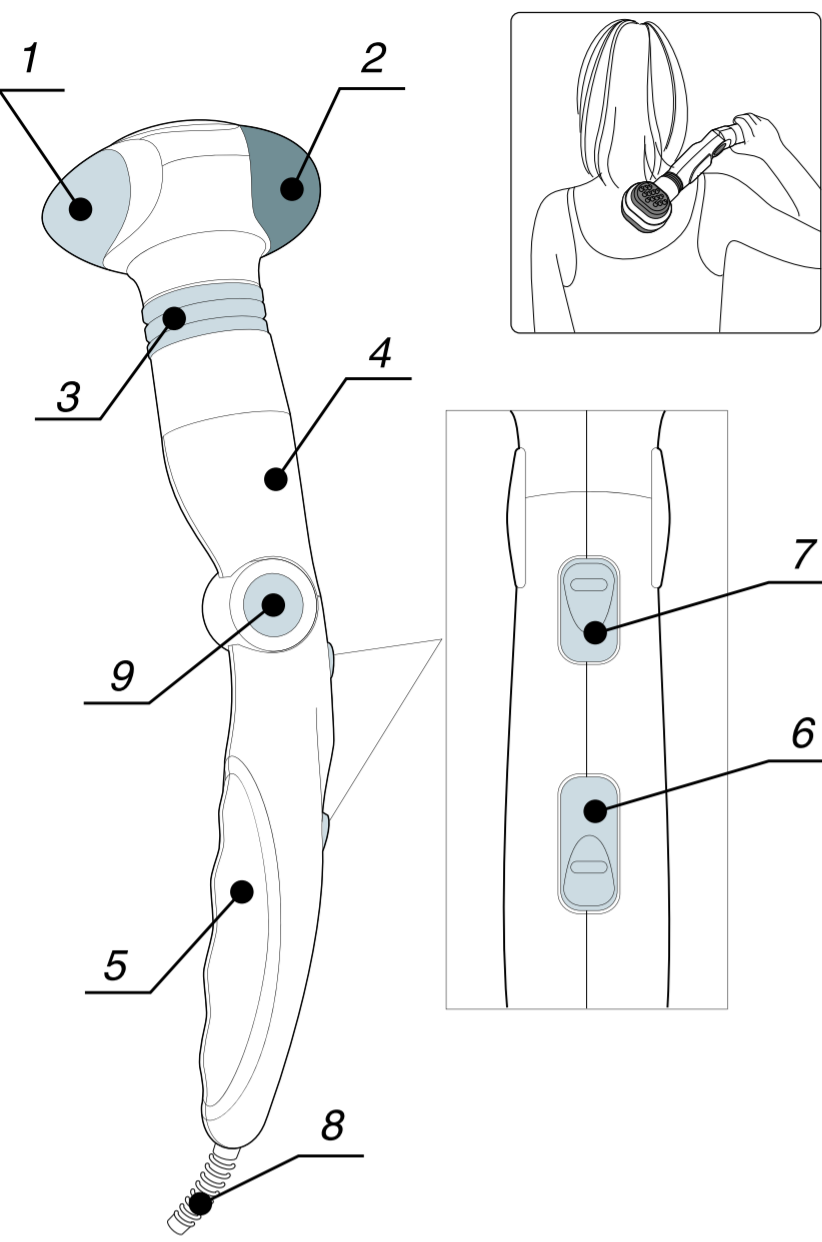**CE** This product conforms to the EMC-Requirements as laid down by the Council Directive 89/336/EEC and to the Low Voltage Regulation (73/23 EEC)

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